

## Culturing Bacteria from *Caenorhabditis elegans* Gut to Assess Colonization Proficiency

Facundo Rodriguez Ayala\*, Sebastián Cogliati, Carlos Bauman, Cecilia Leñini,  
Marco Bartolini, Juan Manuel Villalba, Federico Argañaraz and Roberto Grau

Departamento de Microbiología, Facultad de Ciencias Bioquímicas y Farmacéuticas, Universidad Nacional de Rosario, CONICET, Rosario, Argentina

\*For correspondence: [facundoayala.foncyt@gmail.com](mailto:facundoayala.foncyt@gmail.com)

**[Abstract]** Determining an accurate count of intestinal bacteria from *Caenorhabditis elegans* is one critical way to assess colonization proficiency by a given bacteria. This can be accomplished by culturing appropriate dilutions of worm gut bacteria on selective or differential agarized media. Because of the high concentration of bacteria in gut worm, dilution is necessary before plating onto growth media. Serial dilutions can reduce the concentration of the original intestinal sample to levels low enough for single colonies to be grown on media plates, allowing for the calculation of the initial counts of bacteria in the intestinal sample.

**Keywords:** *Caenorhabditis elegans*, *Bacillus subtilis*, *Escherichia coli*, Culture

**[Background]** Animals rarely live in isolation but rather exist in association with microorganisms. The more characteristic host-microbe interaction in nature is the symbiotic relationship between host and intestinal microbiota (Rosenberg and Zilber-Rosenberg, 2011). In mammals, host-microbe symbiotic interactions mainly occur along mucosal surfaces, with the most important one being the intestinal mucosa. When freshly isolated from the wild, *C. elegans* often harbours a diverse bacterial flora in its gut lumen, reminiscent of the microbial communities of higher organisms (Duveau and Felix, 2012; Bumbarger *et al.*, 2013). By contrast, in the laboratory *C. elegans* is typically maintained in the presence of single bacterial strain (Brenner, 1974). Most often, this is the Gram-negative bacterium *Escherichia coli*. However, other species are sometimes used, such as the Gram-positive *Bacillus subtilis* (Garsin *et al.*, 2003). An adult worm contains approximately 10,000 bacterial cells, a number 10-times greater than that of host worm somatic cells (Portal-Celhay and Blaser, 2012): perhaps coincidentally, this microbiota-to-host cell ratio is similar to that found in humans.

Different strategies are being used to measure the colonization proficiency of *C. elegans* gut by a bacterium such as fluorescein isothiocyanate (FITC)-labelled bacteria, bacteria expressing a fusion reporter (green fluorescent protein [GFP] or  $\beta$ -galactosidase). However, the more accurate method is to measure the number of CFU isolated from the worm intestine. In this protocol, we show how to isolate and count *E. coli* and *B. subtilis* strains from *C. elegans* gut. In the case of *B. subtilis*, we also show how to distinguish vegetative forms from highly resistant spores formed by this bacterium inside the *C. elegans* gut.

## **Materials and Reagents**

1. Pipette tips 2-200  $\mu$ l Eppendorf® epT.I.P.S. (Eppendorf, catalog number: 022492039)
2. Pipette tips 50-1,000  $\mu$ l Eppendorf® epT.I.P.S. (Eppendorf, catalog number: 022492055)
3. Petri dishes 60 x 15 mm 500/cs (Fisher Scientific, catalog number: FB0875713A)
4. Petri dishes 35 x 10 mm 500/cs (Fisher Scientific, catalog number: FB0875711YZ)
5. Corning® 15 ml centrifuge tubes (Corning, catalog number: 430791)
6. Eppendorf® Safe-Lock 1.5 ml microcentrifuge tubes (Eppendorf, catalog number: 022363204)
7. Toothpick
8. 99.95 % Platinum, 0.05 % Iridium wire (3 ft/pk) (Tritech Research, catalog number: PT-9901)
9. OP50 *E. coli* bacteria (University of Minnesota, *C. elegans* Genetics Center, MN)
10. Experimental and control *C. elegans* strains (University of Minnesota, *C. elegans* Genetics Center, MN)
11. Sodium chloride (NaCl) (Sigma-Aldrich, catalog number: S7653)
12. Bacto peptone (BD, Bacto™, catalog number: 211677)
13. Agar (Sigma-Aldrich, catalog number: A1296)
14. Hypochlorite (Sigma-Aldrich, catalog number: 13440)
15. Commercial Bleach 60 g/L (DROGUERÍA INDUSTRIAL SAN JUAN, <http://www.sanjuandrogueria.com>)
16. Triton X-100 (Sigma-Aldrich, catalog number: X100)
17. Luria broth (Sigma-Aldrich, catalog number: L3522)
18. Luria broth with agar (Sigma-Aldrich, catalog number: L2897)
19. Lysozyme from chicken egg white (Sigma-Aldrich, catalog number: L6876)
20. Sodium phosphate dibasic ( $\text{Na}_2\text{HPO}_4$ ) (Sigma-Aldrich, catalog number: S3264)
21. Potassium phosphate monobasic ( $\text{KH}_2\text{PO}_4$ ) (Sigma-Aldrich, catalog number: P5655)
22. Cholesterol (Sigma-Aldrich, catalog number: C8667)
23. 100% ethanol (Sigma-Aldrich, catalog number: E7023)
24. Magnesium sulfate heptahydrate ( $\text{MgSO}_4 \cdot 7\text{H}_2\text{O}$ ) (Sigma-Aldrich, catalog number: M1880)
25. Calcium chloride dihydrate ( $\text{CaCl}_2 \cdot 2\text{H}_2\text{O}$ ) (Sigma-Aldrich, catalog number: C3881)
26. Potassium phosphate dibasic ( $\text{K}_2\text{HPO}_4$ ) (Sigma-Aldrich, catalog number: P2222)
27. Sodium hydroxide (NaOH) (Sigma-Aldrich, catalog number: S8045)
28. Levamisole hydrochloride (Sigma-Aldrich, catalog number: L0380000)
29. Nematode growth medium (NGM) (see Recipes)
30. M9 buffer (see Recipes)
31. 5 mg/ml cholesterol (see Recipes)
32. 1 M  $\text{MgSO}_4$  (see Recipes)
33. 1 M  $\text{CaCl}_2$  (see Recipes)
34. Phosphate buffer (see Recipes)
35. 1 N NaOH (see Recipes)

36. 25 mM Levamisole (see Recipes)

### **Equipment**


1. Erlenmeyer flask (Fisher Scientific, catalog number: FB5006000)
2. Pipetor (Gilson, catalog number: F167300)
3. Worm pick. Worm picks can either be purchased (Genesee Scientific, catalog number: 59-AWP) or made in the lab as described in Wollenberg *et al.*, 2013
4. Pasteur glass pipette (Fisher Scientific, catalog number: 22-378893)
5. Autoclave (Tuttnauerusa, model: 6690)
6. Stirring hotplate (Corning, catalog number: 6795-620)
7. Centrifuge (Eppendorf, model: 5430)
8. Tabletop centrifuge (Eppendorf, model: 5424)
9. Pellet pestle (Kimble Chase Life Science and Research Products, catalog number: 7495211590)
10. Refrigerated incubator (Thermo Fisher Scientific, Thermo Scientific™, model: Heratherm™ General Protocol Microbiological Incubators, catalog number: 51028064)
11. Bunsen burner (Humbolt, catalog number: H-5870)
12. Dissecting stereomicroscope (Olympus, model: SMZ645)
13. Incubators for stable temperature (AQUA® LYTIC incubator 20 °C)
14. Freezers (-20 °C; So-Low Environmental Equipment) (Siemens, model: C85-22)
15. Water bath

### **Procedure**

#### A. Preparation of a synchronized nematode population

1. For NGM plate's preparation, mix 3 g NaCl, 2.5 g Bacto peptone and 17 g agar in Erlenmeyer flask. Add to 1 L of dH<sub>2</sub>O. Autoclave 121 °C for 20 min. Cool flask in 55 °C water bath for 15 min. Add 1 ml 1 M CaCl<sub>2</sub>, 1 ml 5 mg/ml cholesterol in ethanol, 1 ml 1 M MgSO<sub>4</sub> and 25 ml 1 M KPO<sub>4</sub> (Phosphate) buffer (see Recipes). Swirl to mix well. Using sterile procedures dispense the NGM solution into 60 x 15 mm Petri plates. Fill plates 2/3 full of agar.
2. Pick 10 young adults onto an *E. coli* OP50 seeded 60 mm NGM plate (see Recipes). Allow the worms to grow 2-3 days to ensure at least 100 gravid adult worms and adequate number of eggs on the plate.
3. Once reached this worm number, pour 3 ml of M9 buffer (see Recipes) onto the plate and gently swirl it to dislodge the worms. Repeat this procedure one more time using 2 ml of M9 buffer. One plate is enough for bleaching.
4. Using a pipette, transfer the worms to a 15 ml conical tube.
5. Centrifuge for about 1 min at 252 x g to pellet the worms.
6. Aspirate most of the M9 without disturbing the worm pellet.


7. Add about 1 ml of 3% hypochlorite solution and 2.5 ml of 1 N NaOH (see Recipes) to the tube.
  8. Vortex the tube for approximately 5 min or until you see a decrease in the number of intact adult worms. Do not bleach for much longer than this or you will kill the eggs.
  9. Once most of the bodies have dissolved, neutralize the alkaline hypochlorite by adding M9 buffer until 15 ml of the conical tube and centrifuge at 252 x g for 1 min.
  10. Aspirate most of the supernatant solution without disturbing the egg pellet.
  11. Add about 15 ml of M9 to the tube and mix well.
  12. Centrifuge again at 252 x g for 1 min.
  13. Aspirate most of the M9 without disturbing the egg pellet.
  14. Repeat steps A11-A13 at least one more time.
  15. Add about 5 ml of fresh M9 and agitate to resuspend the egg-enriched pellet. Incubate overnight at room temperature with gentle rocker stirrer. The eggs will hatch and the animals will arrest as starved L1 larvae.
  16. Transfer between 100 and 200 starved L1 larvae to a 60 mm NGM plate with *E. coli* OP50 cells ( $1 \times 10^5$  cells per plate) or spores ( $1 \times 10^5$  spores per plate) (see Note 2) of each *B. subtilis* strain and allow the worms to grow 2 days until they reach L4/adult stage.
- B. Preparation of nematode samples for intestinal bacterial counting
1. Use an eyebrow hair and collect 50 adult animals in M9 buffer in a 1.5 ml tube (see Note 1).
  2. Use the same number of animals (worms) fed on each bacterial strain during sample preparation.
  3. Treat the worms with 0.25 mM levamisole to induce temporal paralysis.
  4. Superficially sterilize with 3% commercial Bleach for 5 min (see Note 5).
  5. Wash at least three times with M9 buffer.
  6. After the worms were surface sterilized, worms devoid of outside bacteria are disrupted using a pellet pestle and 1% Triton X-100 (see Note 6).
  7. Centrifuge at 14,800 x g for 10 min, remove the supernatant and resuspend the pellet of each tube in 500  $\mu$ l M9 buffer (see Note 4).
- C. Preparation of sample dilutions (Figure 1)
1. To begin the procedure, take 50  $\mu$ l of the cell suspensions and add to 450  $\mu$ l of deionized water. Shake the suspension well, and label as '10<sup>-1</sup>'.
  2. Before the sample settles, remove 50  $\mu$ l of the cell suspension with a sterile pipette and transfer it to a 450  $\mu$ l-deionized water blank. Vortex the sample, and label as '10<sup>-2</sup>'.
  3. Repeat this dilution step, using 50  $\mu$ l of the previous suspension and a 450  $\mu$ l-deionized water blank. Label these sequentially as tubes '10<sup>-3</sup>'.

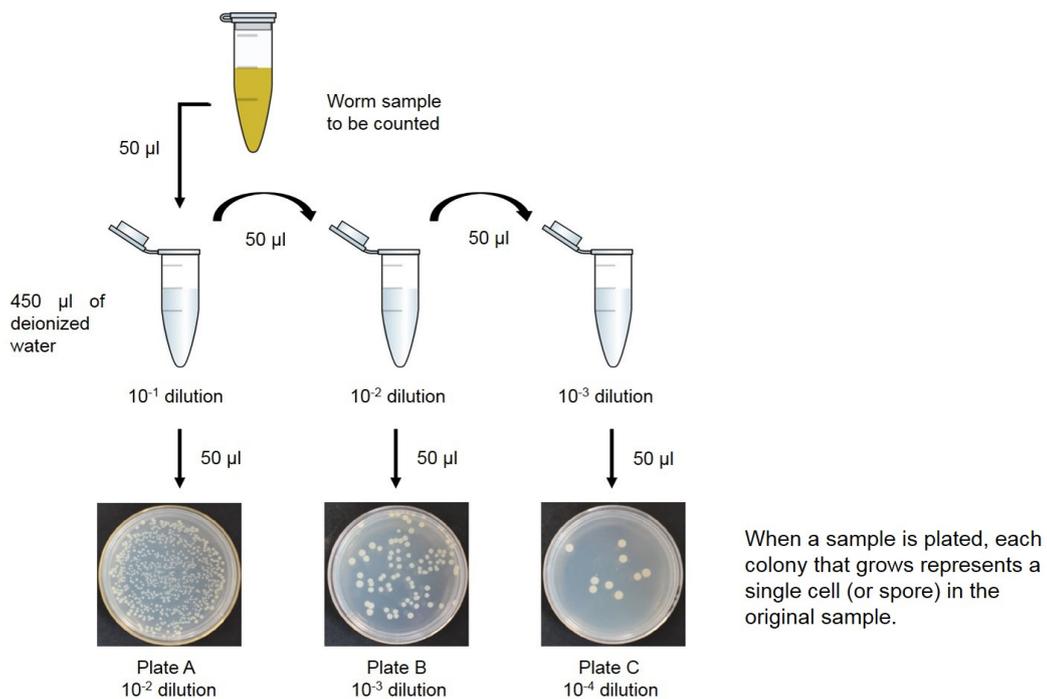

**Figure 1. Diagram of how to prepare bacterial dilutions.** Serial dilution of an initial culture to obtain solutions that are 1/10th, 1/100th and 1/1,000th the concentration of the initial sample ( $10^{-1}$ ,  $10^{-2}$  and  $10^{-3}$  dilutions, respectively). For each dilution, tubes initially have 450 µl of water in them so that the final volume will be 500 µl after the addition of 50 µl of culture.

D. Making spread plates for bacterial culture

- To get bacterial colonies from worm gut isolated *E. coli* and vegetative *B. subtilis* (see Figure 2A), take three LB agar plates (see Note 3) and label them as A, B, and C. Vortex samples  $10^{-1}$ ,  $10^{-2}$ , and  $10^{-3}$ , and pipette 50 µl onto each plate. This increases the dilution value further, by a factor of ten (A =  $10^{-2}$ , B =  $10^{-3}$ , C =  $10^{-4}$ ) (see Figure 1).

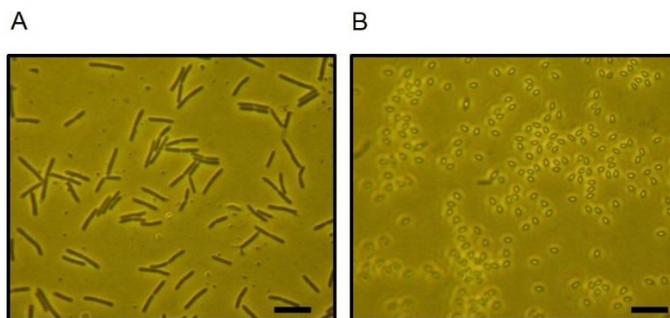

**Figure 2. The micrographs show the vegetative cells (A) and spores (B) of *Bacillus subtilis*.** Bar size = 5 µm.

- Next, dip a glass spreader into ethanol. Place the spreader in a flame for a few seconds to ignite and burn off the ethanol. This will sterilize the spreader.

3. Hold the spreader above the first plate until the flame is extinguished. Open the plate quickly, holding the lid close by. Touch the spreader to the agar away from the inoculum (Inoculum = cells used to begin a culture) to cool, and then spread the drop of inoculum around the surface of the agar until traces of free liquid disappear.
4. Re-flame the spreader and repeat the process with the next plate, working quickly so as not to contaminate the agar with airborne organisms.
5. Incubate the bacteria plates at 37 °C for 16-18 h. Make sure the plates are inverted during the incubation to prevent drops of moisture from condensation from falling onto the agar surface.

#### E. Making spread plates for spore culture

1. Heat-treat the tubes with *B. subtilis* samples labeled as  $10^{-1}$ ,  $10^{-2}$ , and  $10^{-3}$  for 20 min at 80 °C to kill vegetative cells. Vortex spores samples (see Figure 2B) and pipette 50  $\mu$ l onto each plate. This increases the dilution value further, by a factor of ten ( $A = 10^{-2}$ ,  $B = 10^{-3}$ ,  $C = 10^{-4}$ ).
2. Next, dip a glass spreader into ethanol. Place the spreader in a flame for a few seconds to ignite and burn off the ethanol. This will sterilize the spreader.
3. Hold the spreader above the first plate until the flame is extinguished. Open the plate quickly, holding the lid close by. Touch the spreader to the agar away from the inoculum to cool, and then spread the drop of inoculum around the surface of the agar until traces of free liquid disappear.
4. Re-flame the spreader and repeat the process with the next plate, working quickly so as not to contaminate the agar with airborne organisms.
5. Incubate the bacteria plates at 37 °C for 16-18 h. Make sure the plates are inverted during the incubation to prevent drops of moisture from condensation from falling onto the agar surface.

#### F. Bacterial counts

1. After incubation, examine all of the bacteria plates carefully, and note differences in colony size and shape.
2. Count and record the number of bacterial colonies. Only count and use for calculation, plates with 30-300 colonies per plate. A number of colonies per plate below 30 is not statistically significant whereas a number higher than 300 colonies per plate is hard to count, and it is possible a colony arise from more than one cell.
3. The number of CFU/worm is calculated as follows:  
 Number of CFU/worm in plate A = (number of colonies  $\times 10^2$ )/50  
 Number of CFU/worm in plate B = (number of colonies  $\times 10^3$ )/50  
 Number of CFU/worm in plate C = (number of colonies  $\times 10^4$ )/50

## **Data analysis**

1. For each strain or condition, use at least 50-100 animals to obtain results that are more accurate.
2. For each experiment, the same number of nematodes should be examined for each strain and condition.
3. Each assay should be run for triplicate and repeated at least three times (Duveau and Felix, 2012).
4. Use the Student's *t*-test with a significance cut-off level of  $P < 0.05$  for comparisons between two groups.
5. Use the one-factor (ANOVA) variance analysis and correct by the post hoc Bonferroni test for multiple comparisons.

## **Notes**

1. Take a toothpick and glue an eyebrow hair to the tip of it. Let it dry at room temperature. Then, use this tool to pick nematodes. Before using the eyebrow hair always, sterilize it by using ethanol.
2. To obtain pure spores, the heat-treated culture (20 min at 80 °C to kill vegetative cells in a water bath) was treated three times with lysozyme, washed each time with cold water and centrifuged until 100% of the culture consisted of phase-bright spores.
3. For the LB plates used to seed *Bacillus subtilis*, we strongly recommend to dry plates at 45 °C for 20 min or 37 °C for 40 min before to drop the dilutions onto each plate to avoid the sliding movement of this bacterium in solid surface with water drops from condensation.
4. Samples should be used immediately or store for 1 week at -20 °C.
5. To surface sterilize, resuspend 50 worms in 500 µl M9 buffer, add 0.25 mM levamisole (5 µl of 25 mM levamisole stock solution) and 3 % commercial Bleach (stock solution). Mix samples by gentle inversion for 15 min. Then, wash the worms with M9 buffer at least 3 times.
6. After the worms are surface sterilized, pour the 200 µl worm suspension in M9 buffer in a 1.5 ml microcentrifuge tube. Add the extraction solution of choice (1% Triton X-100), and the pestle is used to grind the sample 10 times or until a creamy white to yellow suspension is obtained. After centrifugation, cellular debris and bacteria remain firmly in the bottom of the conical tube, and the supernatant is easily pipetted out.

## **Recipes**

1. Nematode growth medium (NGM)  
Dissolve 3 g NaCl, 2.5 g Bacto peptone and 17 g agar to 1 L of dH<sub>2</sub>O  
Autoclave  
Store at room temperature

2. M9 buffer  
Dissolve 5 g NaCl, 6 g Na<sub>2</sub>HPO<sub>4</sub> and 3 g KH<sub>2</sub>PO<sub>4</sub> to 1 L of dH<sub>2</sub>O  
Autoclave  
Store at room temperature
3. 5 mg/ml cholesterol  
Dissolve 0.25 g of cholesterol in 50 ml of 100% ethanol  
Do not autoclave  
Store at room temperature
4. 1 M MgSO<sub>4</sub>  
Dissolve 6 g MgSO<sub>4</sub> heptahydrate in 50 ml of dH<sub>2</sub>O  
Autoclave  
Store at room temperature
5. 1 M CaCl<sub>2</sub>  
Dissolve 5.55 g CaCl<sub>2</sub> dihydrate in 50 ml of dH<sub>2</sub>O  
Autoclave  
Store at room temperature
6. Phosphate buffer  
Dissolve 10.7 g K<sub>2</sub>HPO<sub>4</sub> and 32.5 g KH<sub>2</sub>PO<sub>4</sub> to 300 ml of dH<sub>2</sub>O  
Adjust pH to 6.0  
Autoclave  
Store at room temperature
7. 1 N NaOH  
Dissolve 2 g NaOH in 50 ml of dH<sub>2</sub>O
8. 25 mM levamisole  
Dissolve 0.3 g levamisole in 50 ml of dH<sub>2</sub>O

### **Acknowledgments**

This work was supported by CONICET (Consejo Nacional de Investigaciones Científicas y Técnicas) and FONCyT (Fondo para la Investigación Científica y Tecnológica) with the aid of the Pew Latin-American Program in Biological Sciences (Philadelphia, PA, USA), the Fulbright Committee (Washington, DC, USA) and former Fundación Antorchas (Buenos Aires, Argentina). We modified the media and NGM plate's preparation from Stiernagle T. Maintenance of *C. elegans*. WormBook. 2006 11:1-11. We modified *C. elegans* synchronization from Montserrat Porta-de-la-Riva *et al.* (2012). Basic *Caenorhabditis elegans* Methods: Synchronization and Observation. Journal of Visualized Experiments; 64, e4019, 1-9.

## References

1. Brenner, S. (1974). [The genetics of \*Caenorhabditis elegans\*](#). *Genetics* 77(1): 71-94.
2. Bumbarger, D. J., Riebesell, M., Rodelsperger, C. and Sommer, R. J. (2013). [System-wide rewiring underlies behavioral differences in predatory and bacterial-feeding nematodes](#). *Cell* 152(1-2): 109-119.
3. Duveau, F. and Felix, M. A. (2012). [Role of pleiotropy in the evolution of a cryptic developmental variation in \*Caenorhabditis elegans\*](#). *PLoS Biol* 10(1): e1001230.
4. Garsin, D. A., Villanueva, J. M., Begun, J., Kim, D. H., Sifri, C. D., Calderwood, S. B., Ruvkun, G. and Ausubel, F. M. (2003). [Long-lived \*C. elegans daf-2\* mutants are resistant to bacterial pathogens](#). *Science* 300(5627): 1921.
5. Porta-de-la-Riva, M., Fontrodona, L., Villanueva, A. and Cerón, J. (2012). [Basic \*Caenorhabditis elegans\* methods: synchronization and observation](#). *J Vis Exp* 64: e4019.
6. Portal-Celhay, C. and Blaser, M. J. (2012). [Competition and resilience between founder and introduced bacteria in the \*Caenorhabditis elegans\* gut](#). *Infect Immun* 80(3): 1288-1299.
7. Rosenberg, E. and Zilber-Rosenberg, I. (2011). [Symbiosis and development: the hologenome concept](#). *Birth Defects Res C Embryo Today* 93(1): 56-66.
8. Stiernagle, T. (2006). [Maintenance of \*C. elegans\*](#). *WormBook* 11:1-11.
9. Wollenberg, A. C., Visvikis, O., Alves, A. F. and Irazoqui, J. E. (2013). [Staphylococcus aureus killing assay of \*Caenorhabditis elegans\*](#). *Bio-protocol* 3(19): e916.